Summary


Summary	1
Description	2
English	3
Français	5
Español	7
Deutsch	9
Italiano	11
Português	13
Čeština	15
Nederlands	17
Polski	19
Glossary	21

DENTAL MONITORING


by DentalMind

www.dental-monitoring.com

1

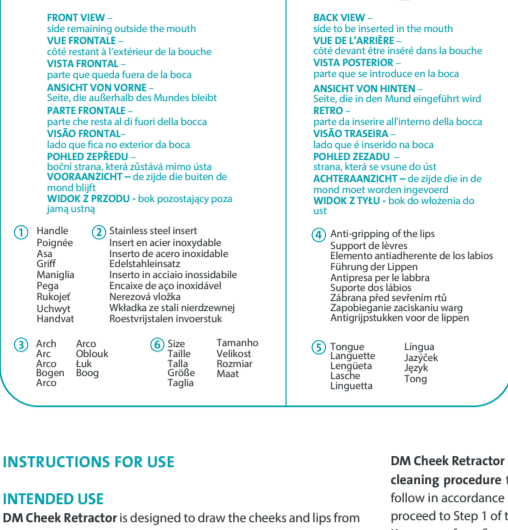DESCRIPTION / DESCRIPCIÓN / DESCRIZIONE / DESCRICÃO / POPIS / BESCHRIJVING / OPIS

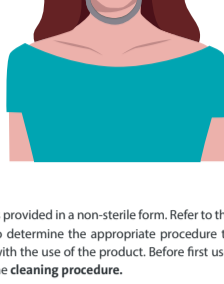2

INSTRUCTIONS FOR USE

INTENDED USE

DM Cheek Retractor is designed to draw the cheeks and lips from the buccal and labial surfaces of the teeth and gums during scanning procedures performed with products marketed by Dental Monitoring. It is indicated for use for patients over the age of 6 and is to be used in:

- professional healthcare facilities, such as a dental practice, or;
- a non-healthcare environment, such as the patient's own home.

DM Cheek Retractor is available in three different sizes to accommodate different morphologies:

- Size S (REF 002B03004)
- Size M (REF 002B02004)
- Size L (REF 002B01004)

Appropriate sizing should be determined within a healthcare facility.

DM Cheek Retractor is supplied in a transparent sleeve.

CONTRAINDICATIONS

DM Cheek Retractor is not intended for surgical use, nor to take intra-oral pictures with products other than those marketed by Dental Monitoring.

DM Cheek Retractor should not be used in the event of injured skin or injured mucous membrane.

DM Cheek Retractor should not be used by patients under the age of 6.

WARNINGS, PRECAUTIONS

en

pursue use.

Note: If exclusively used by one patient, do not proceed to Step 2 and Step 3 and place the product back in its sleeve. If used by multiple patients, please proceed to Step 2 and Step 3.

Step 2: Disinfection

- Disinfect with an acid-based biocide such as Anioxyde® 1000, Bioxal® M, Dynacide® PA, Nu Cidex®.

- Fully submerge the retractor in the solution, and cover as prescribed, please allow the dilution and soaking time as advocated by the manufacturer.

- Remove from the solution and check if there are any noticeable debris. Please gently brush with a soft bristle brush if needed.

- Rinse heavily with water to clear away the chemicals.

- Dry with a clean dry and soft, lint-free cloth or with sterilized cleaning cloths or with previously filtered compressed air.

Step 3: Autoclave Sterilization

- Insert the dried and cleaned retractor in a sterilization pouch.

- Sterilize with steam water at 134°C (273°F) / 2.1 bars / 3 min.

- Store in a dry and clean place without removing the retractor from the sterilization pouch.

DM Cheek Retractor can be autoclaved up to 20 times.

DM Cheek Retractor should be sterilized individually, meaning by inserting only one retractor per sterilization pouch. Use appropriate sizing with regard to the sterilization pouch; the retractor should not be compressed in any way once in the

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it

DM Cheek Retractor is provided in a non-sterile form. Refer to the **cleaning procedure** to determine the appropriate procedure to follow in accordance with the use of the product. Before first use, proceed to Step 1 of the **cleaning procedure**.

Keep away from fire.

Do not clean in a dishwasher.

Do not clean in a dry-heat sterilizer

PUTTING ON YOUR RETRACTOR

In this section, the numbers refer to the device parts as presented in the **description**.

- Hold the retractor using both handles ①, arch down ③ with the stainless-steel inserts away from you ②.

- Slightly fold the retractor towards you by bringing closer the left and right side.

- Insert one side first, then the other. Make sure the tongues are correctly inserted into the mouth ⑤.

- Check if the retractor is correctly applied: teeth and gums should be visible without experiencing pain or discomfort.

In the event of pain, refer to your healthcare facility to check you have the appropriate retractor size.

CLEANING PROCEDURE

Step 1: Basic cleaning

- Clean with soap and cold water.

- Gently brush the surface and make sure it is completely clean.

- Rinse with cold water.

- Dry with a clean dry and soft, lint-free cloth or with previously filtered compressed air.

If you notice signs of over-wear or damage such as cracks, do not

autoclave.

STORAGE

The retractor should be stored in its sleeve or in the sterilization pouch away from solar radiation and moisture.

DISPOSAL

Once your retractor is no longer required, dispose of it in accordance with local regulations.

LIST OF REQUIRED NON-PROVIDED PRODUCTS

- Dental Monitoring ScanBox (REF 001A08001): optional accessory

- Disinfectant with an acid-based biocide (if used by multiple patients)

PERFORMANCE STUDIES

DM Cheek Retractor biocompatibility was assessed and the product was demonstrated to be non-cytotoxic, non-sensitizing and non-irritating.

WARRANTY

One year from the date of original purchase.

TRAVEL

In the event of travel, the retractor should be kept in its sleeve. Please make sure the retractor is not too tightly packed to prevent any deformation during transport.

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